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DEAR CUSTOMER,

Congratulations on your new hearing systems. You have chosen sophisticated, high-quality products. By purchasing these hearing systems, you have taken an important step to enable yourself to hear well again. We are sure these hearing system will soon become indispensable to you and make your life safer and more enjoyable. This user manual will help you to familiarize yourself with your new hearing systems.

If you require further information, please contact your hearing system specialist, who will be happy to help.


THE HEARING SYSTEMS	24-25
HANDLING THE BATTERIES	26
– Battery type	26
– Inserting the battery	26
– Battery warning signal	26
SWITCHING THE HEARING SYSTEM ON AND OFF	27
– Audiomatic (P) power-on delay	27
PROGRAM SELECTION	28
VOLUME CONTROL	29
INSERTING THE HEARING SYSTEM	30
PROTECTION/FILTER	31-32
– Microphone filter	31
– Cerumen protection filter (wax guard)	31-32
WIRELESS	33
– Functionality	33
– Technical data for wireless hearing systems	33
GENERAL SAFETY NOTES	34-38
DISPOSAL/CARE/CE-MARK	39
– Disposal instructions	39
– Cleaning the hearing system ¹	39
– CE-mark	39
CONFIGURATION OF YOUR HEARING SYSTEMS	40-41

THE HEARING SYSTEMS


Icon


sina)xs

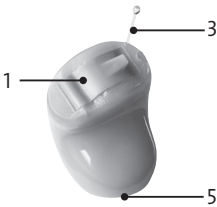
Sina

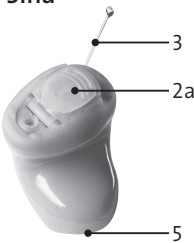Vega

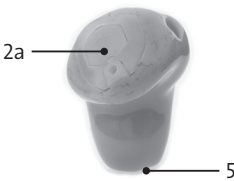Aumea

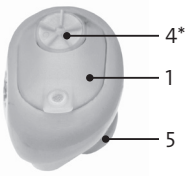Please tick customized hearing system.

1	Battery compartment lid
2	Battery compartment lid with integrated push button*, can be programmed by your hearing system specialist as <ul style="list-style-type: none">- Program button- Program button with on/off function * Depends on the construction of the unit
2a	Battery compartment lid with integrated push button, can be programmed by your hearing system specialist as <ul style="list-style-type: none">- Program button- Program button with on/off function- Volume control (only applies to hearing systems with wireless functionality)
3	Removal line for removing the hearing system from the ear
4	Volume control; * optionally also with program button
5	Wax guard (red/blue)

HANDLING THE BATTERIES

BATTERY TYPE

Please read the information under “Hearing systems passport” on the back page.

When using zinc-air batteries, remove the protective film from the plus side before inserting. The zinc-air battery then needs a few minutes to reach full power.

INSERTING THE BATTERY

Place your fingernail in the small groove or under the small lug on the edge of the battery compartment and open the lid. If there is a spent battery inside, remove it. Insert the new battery (see figure) and close the battery compartment.

The battery compartment is designed to prevent the battery from being inserted the wrong way round. Please do not use force when closing it.

BATTERY WARNING SIGNAL

When battery power is low, the hearing system becomes gradually quieter at first. Then a battery warning signal is heard. You must then insert a new battery. The hearing system will continue to work for approx. 10 minutes after the battery warning signal has sounded.

SWITCHING THE HEARING SYSTEM ON AND OFF

The battery compartment also serves as an on/off switch. Your hearing system is switched on when the battery compartment is fully closed. The hearing system is switched off when the battery compartment is opened again.

If your hearing system has a push button integrated into the battery compartment lid, this can also be used to switch it on and off. Ask your hearing system specialist for details, and refer to the overview on page 41. Some hearing systems can optionally be switched on and off with a remote control. Your hearing system specialist will be happy to answer any questions you may have on this. Always switch off your hearing system after use. If you are not going to use your hearing system for a long time, please take the battery out.

AUDIOMATIC (P) POWER-ON DELAY

If your hearing system features Audiomatic automatic power-on delay, this can be set by your hearing system specialist. After switching your hearing system on, you have up to 12 seconds to insert the system into your ear without feedback whistling.

(P) = Patent Registered

PROGRAM SELECTION

If your hearing system features a push button integrated into the battery compartment lid, you can select from up to 6 different hearing programs, depending on the model.

You can switch between the different hearing programs simply by pressing the program button. You will also hear signal tones that provide acoustic confirmation of the switch to a new program. Ask your hearing system specialist for details and refer to the overview on page 41. With selected hearing systems, the program can also be changed by means of optional remote controls. Your hearing system specialist will be happy to provide you with further information.

VOLUME CONTROL

BY MEANS OF A VOLUME CONTROL KNOB (AUMEA)

If your hearing system has a volume control knob, use this to adjust the volume. Increase the volume by lightly pressing the volume control knob and turning it forwards with the tip of your index finger. Turn the knob backwards to turn the volume down.

WITH THE PUSH BUTTON (SINA/VEGA WITH WIRELESS FUNCTIONALITY)

Your hearing system specialist can program hearing systems with wireless functionality so that the volume control function is operated with the push button. In this case, proceed as follows: Press the push button of one hearing system briefly to increase the volume. Press the push button of the other hearing system briefly to reduce the volume. To see how the volume control is distributed over both sides, refer to the overview on p. 41.

AUTOMATICALLY

With these systems, the amplifier automatically adjusts to the hearing situation and regulates the volume. This protects the ear from volumes that are too loud. The volume cannot be adjusted manually. When fitting your hearing system, your hearing system specialist will have already set it so that it corresponds to your personal hearing perception.

With selected hearing systems, it is also possible to adjust the volume by means of optional remote control units. Please direct any questions you may have on this option to your hearing system specialist.

INSERTING THE HEARING SYSTEM

Insert your hearing system into your ear so that the tip with the sound opening points in the direction of the auditory canal.

The color of the type designation (laser engraving) or the color of the wax guard indicates which is the right hearing system and which the left:

blue = left red = right

If you have an Icon hearing system, you will be able to distinguish the hearing systems for the right and left ears by the white marks on the hearing system housing:

R = right L = left

This mark must face upwards when the hearing system is inserted.

Pull your earlobe a little with your free hand, and at the same time push your hearing system into your ear with the other hand. The housing, which has been adapted to fit your auditory canal, slides into the correct position. When removing the hearing system, perform the sequence in reverse. If your hearing system has a removal line, the hearing system can simply be removed from the ear using this. Your hearing system specialist will be happy to demonstrate to you how to use your hearing system, if you wish.

PROTECTION/FILTER

MICROPHONE FILTER

The microphone filter provides even more effective protection from soiling and ensures optimum long-term transmission performance.

If your hearing system has a microphone filter, this can be changed when necessary by your hearing system specialist.

CERUMEN PROTECTION FILTER (WAX GUARD)

A cerumen protection filter prevents moisture and cerumen (ear-wax) from penetrating the hearing system. If your hearing system has this filter, it should be replaced at least every two weeks. Please use the red filter for the right hearing system (red lettering) and the blue filter for the left hearing system (blue lettering).


If you experience problems when changing the filter, your hearing system specialist will be happy to help.

1. Remove used wax guard:


with thread

2. Dispose of used wax guard:

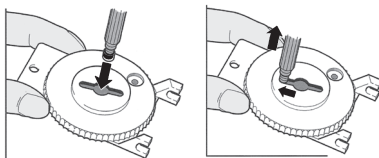
3. Remove new wax guard:

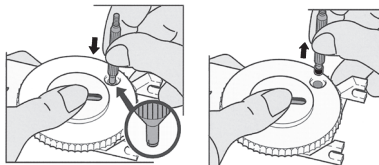without thread

4. Insert new wax guard into the hearing system:

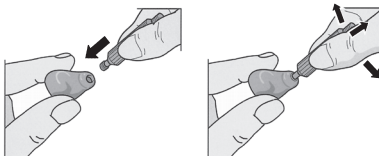WIRELESS

FUNCTIONALITY

If your hearing system has wireless functionality, this supports the following functions:

- Use of a remote control.

- Automatic adjustment of both hearing systems at the same time when you change the system settings, e.g. volume or hearing program.

i The performance of the wireless function can be negatively affected by electromagnetic radiation, for example from computer monitors.

- If disturbance occurs, move further away from the source of the interference.

TECHNICAL DATA FOR WIRELESS HEARING SYSTEMS

Frequency range	3.28 MHz
Pulse duration	325 ms
Transmitting power	42 μ W

GENERAL SAFETY NOTES

i Only use the hearing systems and accessories as described in the user manual.

CAUTION

Please read this user manual carefully and in full. Follow the safety instructions in this document to prevent damage or injury.

WARNING

Possible damage to remaining hearing ability.

- Only wear hearing systems that have been specifically fitted for you.

WARNING

Risk of injury!

- Never use devices with obvious signs of damage. Please return these to the place of sale.

WARNING

Risk of choking if small parts are swallowed.

- Keep hearing systems, batteries, and accessories out of the reach of children.
- If parts are swallowed, go to a doctor or hospital immediately.
- This hearing system is not suitable for children under 36 months or people with mental disabilities. Your hearing system specialist can advise you on selecting a suitable hearing system.

WARNING

Risk of explosion!

- Do not use your hearing systems in environments where there is a risk of explosion (e.g. mines).

WARNING

Electrical interference!

- In environments where the use of electronic or wireless devices is restricted, please check whether your hearing system needs to be switched off.

NOTE

Your hearing systems have been developed in accordance with international electromagnetic compatibility standards. However, they may still encounter interference from electrical devices in the vicinity. If this occurs, move further away from the source of the interference.

NOTE

Leaking batteries will damage the hearing systems.

- Switch off the hearing systems when not in use to conserve battery power.
- Remove batteries if you are not going to use the hearing systems for a long time.

NOTE

Batteries contain hazardous substances that can pollute the environment.

- Please do not dispose of your spent batteries in household waste.
- Please dispose of batteries in accordance with current regulations, or hand them in to your hearing system specialist.

i

Some countries impose certain restrictions on the use of radio devices.

- More information can be obtained from the responsible local authorities.

NOTE

Your hearing systems are sensitive to extreme heat, high humidity, strong magnetic fields (> 0.1 T), X-rays, and mechanical stress.

- Never expose your hearing systems to extreme temperatures or high humidity.
- Do not leave them lying in the sun.
- Do not wear your hearing systems in the shower or when you are applying makeup, perfume, aftershave, hairspray, or sunscreen.
- Do not put your hearing systems in the microwave.

NOTE

Various kinds of strong radiation can damage your hearing systems, such as medical X-rays or magnetic resonance imaging of the head.

Do not wear your hearing systems during these or similar examinations and procedures. Weaker radiation such as that emitted by radio equipment or during security checks at airports does not damage your hearing systems.

DISPOSAL INSTRUCTIONS

Within the European Union, the accessories listed are covered by the European Parliament and Council Directive 2002/96/EC from 27 January 2003 on waste electrical and electronic equipment.

Most recently modified by Directive 2003/108/EC (WEEE).

CLEANING THE HEARING SYSTEM

Various care products are available for cleaning and drying your hearing system. Please ask your hearing system specialist which are most suitable for your hearing system.

Do not use solvents!

Dispose of hearing systems, batteries, and accessories in accordance with the applicable regulations.

CE-MARK

With the CE-mark, Audio Service confirms compliance with European Directive 93/42/EEC on medical devices.

For products with Wireless 2.0, Audio Service confirms compliance with European Directive 99/5/EC (R&TTE) concerning radio equipment and telecommunications terminal equipment.

Subject to technical modifications and errors.

CONFIGURATION OF YOUR HEARING SYSTEMS

Hearing system		
left		right
	Model	
	Battery size	
<input type="checkbox"/> yes <input type="checkbox"/> no	Power-on delay	<input type="checkbox"/> yes <input type="checkbox"/> no
<input type="checkbox"/> yes <input type="checkbox"/> no	Program button	<input type="checkbox"/> yes <input type="checkbox"/> no
Signal tones		
<input type="checkbox"/> on <input type="checkbox"/> off	On/off switch	<input type="checkbox"/> on <input type="checkbox"/> off
<input type="checkbox"/> on <input type="checkbox"/> off	Program change	<input type="checkbox"/> on <input type="checkbox"/> off
<input type="checkbox"/> on <input type="checkbox"/> off	Battery warning signal	<input type="checkbox"/> on <input type="checkbox"/> off
<input type="checkbox"/> on <input type="checkbox"/> off	Volume limit	<input type="checkbox"/> on <input type="checkbox"/> off
<input type="checkbox"/> on <input type="checkbox"/> off	Volume change	<input type="checkbox"/> on <input type="checkbox"/> off
Program	Description of the hearing situation	
1		
2		
3		
4		
5		
6		
	Remote control	
Model		

Hearing system				
left			right	
Push button				
<input type="checkbox"/>	Quick press	Program change	Quick press	<input type="checkbox"/>
<input type="checkbox"/>	Program change with on/off function			<input type="checkbox"/>
	Quick press	Program change	Quick press	
	Hold down for approx .2 sec.	Switch on and off	Hold down for approx .2 sec.	
<input type="checkbox"/>	Quick press	Volume control	Quick press	<input type="checkbox"/>

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Audio Service GmbH · Zeppelinstraße 9 · D-32051 Herford